


# TOC 1000

## User Manual

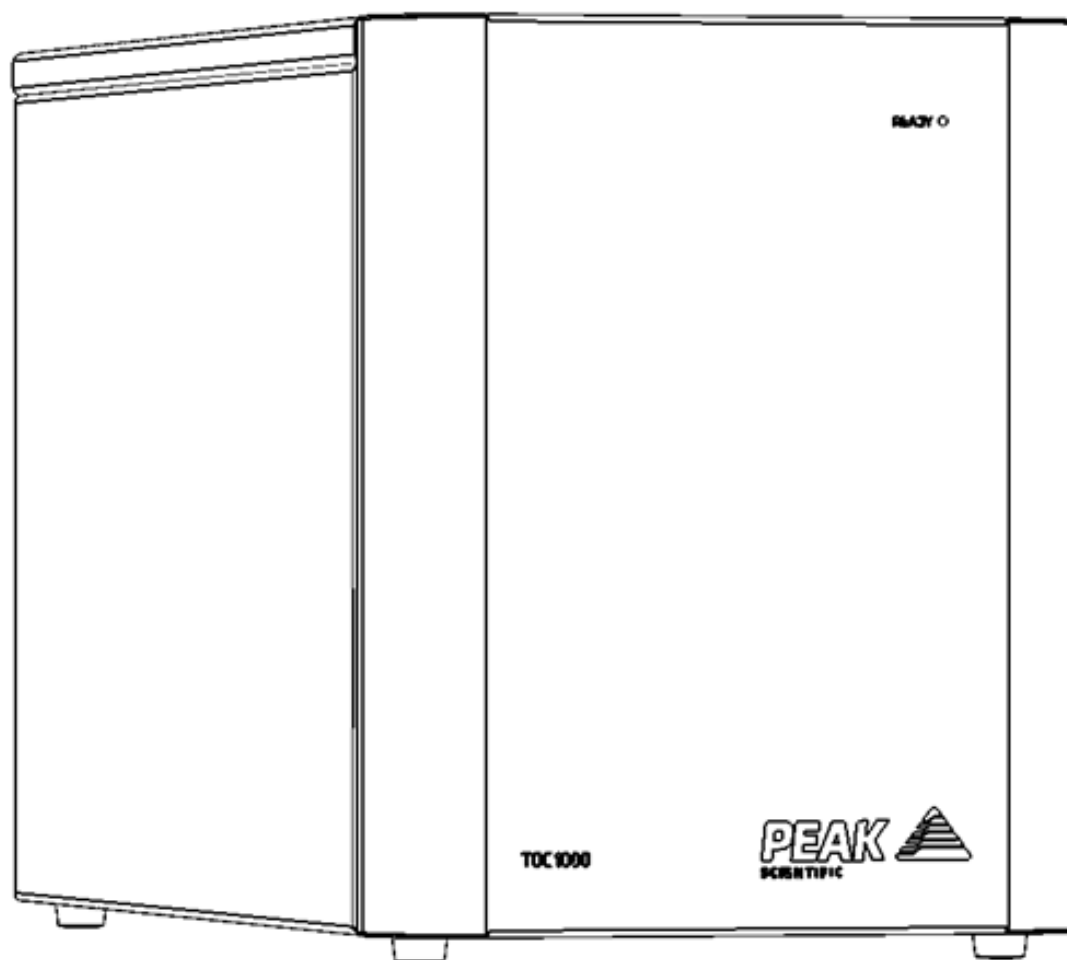

# Contents

Change History	3
How to use this Manual	3
Warranty Statement	4
Safety Notices	5
Symbols	5
Safety Notice to Users	5
EU Declaration of Conformity	6
UK Declaration of Conformity	7
WEEE Compliance Statement	8
European Union (EU) and United Kingdom (UK) Class A Compliance statement	9
FCC Class A Compliance Statement	9
Industry Canada Class A emission compliance statement	9
Korea Communications Commission (KCC) statement	9
CSA Compliance Statement	10
Technical Specification	11
TOC 1000	11
Unpacking	12
Fittings Kit Contents	13
Installation	14
Generator Environment	14
Generator Overview	15
General Dimensions	15
Rear Connections	16
Electrical Connection	17
Air Purity	18
Class 1 Particulate	18
Class 4 Water	18
Class 1 Oil	18
Start-Up Sequence	19
Connecting to the application	20
Tubing Lengths	20
Service Requirements	21
Service Schedule	21
Peak Protected	22
Cleaning	23
Troubleshooting	24

## Change History

Rev	Comment	Name	Date
1	Initial Release	Liam Couttie	18/02/2020
2	Declarations Update	Cleo Denholm	31/08/2021
3	Compressed Air Supply Update	L. Couttie	28/07/2023
4	Spec. Update	L. Couttie	22/03/2024
5	Sticker added on	Cleo Denholm	07/06/2024
6	Compressed Air Note Added	L. Couttie	09/08/2024
7	Spec. Update	L. Couttie	12/11/2024
8	Declarations Update	L. Couttie	26/03/2025

## How to use this Manual

This manual is intended for end users and has been written as a reference document where you can skip to the relevant information.

Users can refer to the contents page to find the relevant information.

Please review each of the following sections carefully.

Thank you for selecting Peak Scientific to meet your gas generation needs, and should you require any further assistance or support please do not hesitate to contact Peak Scientific or the Peak Partner from which you purchased your generator.

## **Warranties and Liabilities**




**Visit:** [www.peakscientific.com/warranty-statement/](http://www.peakscientific.com/warranty-statement/)

# Safety Notices

Peak Scientific Instruments cannot anticipate every possible circumstance which may represent a potential hazard. The warnings detailed within this manual refer to the most likely potential hazards, but by definition cannot be all inclusive. If the user employs an operating procedure, item of equipment or a method of working which is not specifically recommended by Peak Scientific, the user must ensure that the equipment will not be damaged or become hazardous to persons or property.

## Symbols

This manual uses the following symbols to highlight specific areas important to the safe and proper use of the generator.

 WARNING	A WARNING notice denotes a hazard. It calls attention to an operating procedure, process or similar, which if not correctly performed or adhered to, could cause personal injury or in the worst case death. Do not proceed beyond a WARNING notice until the indicated conditions are fully understood or met.
 CAUTION	A CAUTION notice denotes a hazard. It calls attention to an operating procedure, process or similar, which if not correctly performed or adhered to, could cause damage to the generator or the application. Do not proceed beyond a CAUTION notice until the indicated conditions are fully understood or met.
	Caution, risk of electric shock. Ensure power to the generator has been removed before proceeding.

## Safety Notice to Users



These instructions must be read thoroughly and understood before installation and operation of your Peak Generator. Use of the generator in a manner not specified by Peak Scientific MAY impair the SAFETY provided by the equipment.



When handling, operating or carrying out any maintenance, personnel must employ safe engineering practices and observe all relevant local health and safety requirements and regulations. The attention of UK users is drawn to the Health and Safety at Work Act 1974, and the Institute of Electrical Engineers regulations.



If the equipment is used in a manner not specified by the manufacturer, the protection provided by the equipment maybe impaired.

# EU Declaration of Conformity

We Peak Scientific Instruments Ltd.  
Of Fountain Crescent, Inchinnan, Renfrewshire, PA4 9RE

Hereby declare that, this declaration of conformity is issued under the sole responsibility of the manufacturer.

Equipment Type: TOC Air Generator  
Model Designator: TOC 1000

To which this declaration relates, is in conformity with the following applicable EU Directives, harmonized standards, and other normative requirements.

- **Low Voltage Directive 2014/35/EU**  
EN 61010-1: 2010+ A1:2019 Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use.
- **Electromagnetic Compatibility Directive 2014/30/EU**  
EN 61326-1: 2013 Electrical Equipment for Measurement, Control and Laboratory Use – EMC Requirements. (Class A)
- **Restriction on the use of certain hazardous substances in electronic equipment (RoHS) Directive 2011/65/EU as amended by EU 2015/863.**

Signed for and on behalf of Peak Scientific by

Signed: 

Name: Fraser Dunn

Position: Design Engineering Manager  
Peak Scientific Instruments Ltd,  
Inchinnan, Renfrew, Scotland, PA4 9RE, UK.

Date: 14th July 2024



# UK Declaration of Conformity

We Peak Scientific Instruments Ltd.  
Of Fountain Crescent, Inchinnan, Renfrewshire, PA4 9RE

Hereby declare that, this declaration of conformity is issued under the sole responsibility of the manufacturer.

Equipment Type: TOC Air Generator  
Model Designator: TOC 1000

To which this declaration relates, is in conformity with the following applicable UK Statutory Instruments, Standards and other normative requirements.

- **The Electrical Equipment (Safety) Regulations 2016 (SI 2016 / 1101) as amended.**  
BS61010-1:2010+ A1:2019 Safety Requirements for Electrical Equipment for Measurement Control and Laboratory Use.
- **The Electromagnetic Compatibility Regulations 2016 (SI 2016 / 1091) as amended.**  
BS61326-1:2013 Electrical Equipment for Measurement , Control and Laboratory Use  
- EMC Requirements.
- **The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012 (SI 2012 / 3032) as amended.**

Signed for and on behalf of Peak Scientific by

Signed: 

Name: Fraser Dunn

Position: Design Engineering Manager  
Peak Scientific Instruments Ltd,  
Inchinnan, Renfrew, Scotland, PA4 9RE, UK.


Date: 14th July 2024



# WEEE Compliance Statement

The Waste Electrical and Electronic Equipment (WEEE) Regulations SI 2013 No 3113 and or the Waste Electrical and Electronic Equipment (WEEE) Directive 2012/19/EU apply to all electrical and electronic equipment placed on the market in the UK and EU covered by the scope of regulations which can be found in the Government Guidance Notes (PDF) produced by the Department for Business Innovation and skills for the UK and here for Europe.

All PEAK products that are subject to the WEEE directive are compliant with the WEEE marking requirement. Such products are marked with the “crossed-out wheellie bin” symbol (shown below) in accordance with European standard EN50419. All old electrical equipment can be recycled. Please do not dispose of any electrical equipment (including those marked with this symbol) in general rubbish bins. Please contact your dealer or distributor for clarity.

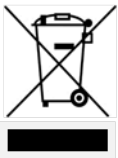



# **EMC Class A Compliance Statements**

## **European Union (EU) and United Kingdom (UK) Class A Compliance statement**

This is a Class A product. In a domestic environment this product may cause radio interference in which case the user may be required to take adequate measures.

## **FCC Class A Compliance Statement**

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.



## **Industry Canada Class A emission compliance statement**

This ISM device complies with Canadian ICES-001 (A).

Cet appareil ISM est conforme à la norme NMB-001 (A) du Canada.

## **Korea Communications Commission (KCC) statement**

이 기기는 업무용(A급)으로 전자파적합기기로서 판매자 또는 사용자는 이 점을 주의하시기 바라며, 가정외의 지역에서 사용하는 것을 목적으로 합니다.

(This is electromagnetic wave compatibility equipment for business (Type A). Sellers and users need to pay attention to it. This is for any areas other than home.)

# CSA Compliance Statement

CSA Group (Canadian Standards Authority) is a Nationally Recognised Testing Laboratory (NRTL), headquartered in Toronto Canada.

They are authorised to evaluate product to both their own and Underwriters Laboratory (UL) standards and certify the product to be in compliance to the relevant standards.

Peak products are certified to the current in force revision of the following standards in order to cover both Canadian and United States requirements for “Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory use, Part 1: general Requirements”.

Canada: CAN/CSA C22.2 No 61010-1-12

United States: UL 61010-1

As a result the products covered by this statement are certified and listed by CSA accordingly and are entitled to carry the CSA mark with both Canadian and United States subscripts , as shown below on the product rating label.



# Technical Specification

## TOC 1000

### Environment

Minimum Operating Ambient Temperature	5°C (41°F)
Maximum Operating Ambient Temperature	35°C (95°F)
Maximum Altitude	2000m
Maximum Relative Humidity	70% @ 35°C
Minimum Storage Temperature*	5°C (41°F)
Maximum Storage Temperature*	35°C (95°F)
Maximum Storage Temperature*	90% non-condensing

### Inlet Conditions

Min/Max Inlet Air Pressure	100psi / 6.9bar - 120psi / 8.3bar
Minimum Air Inlet Flow	6 LPM
Minimum Air Inlet Quality	ISO 8573 - 1:2010 Class 1.4.1†
Maximum Inlet CH <sub>4</sub> Concentration	100ppm CH <sub>4</sub>

### Generator Outlets

Maximum Gas Output Pressure	50 psi / 3.4 bar
Maximum Outlet Flow Rate**	1000 cc/min
Gas Purity	<1ppm CO <sub>2</sub> <1ppm CO <1ppm SO <sub>x</sub> <1ppm NO <sub>x</sub> <0.05ppm CH <sub>4</sub>
Dewpoint (°C)	-60°C / -76°F
Particles	<0.01 µm
Phthalates	Phthalate & BHT Free
Suspended Liquids	None
Gas Outlets	1 x 1/4" BSP x 6mm push-fit connection
Start-Up Time	120 minutes

### Electrical Requirements

Voltage	120V±5% / 230V ±10%
Frequency	120V/ 60Hz / 230V 50/60Hz
Current	120V 4A / 230V 1A
Power Cord (Supplied)	C13 Plug
Circuit Breakers	10A
Pollution Degree	2
Insulation Category	Class 1
Transient Over Voltages	Over Voltage Category II

### General

Dimensions mm (inches) HxWxD	343 x 270 x 530 mm (13.5 x 10.7 x 20.9")
Shipping Dimensions mm (inches) WxDxH	435 x 650 x 400mm (17.1 x 25.5 x 15.7")
Generator Weight Kg (lbs)	27kg (59.6lbs)
Shipping Weight Kg (lbs)	38kg (83.8lbs)
Noise Level	Virtually Silent


\* Note: Please ensure generator is situated in a well ventilated environment.

\*\* Note: Flows in LPM are expressed as normalised volumes at 101.3kPa, 20°C.

# Unpacking

Although Peak Scientific takes every precaution with safe transit and packaging, it is advisable to fully inspect the unit for any sign of transit damage.

Check 'SHOCKWATCH' and 'TIP-N-TELL' labels for signs of rough handling prior to unpacking.

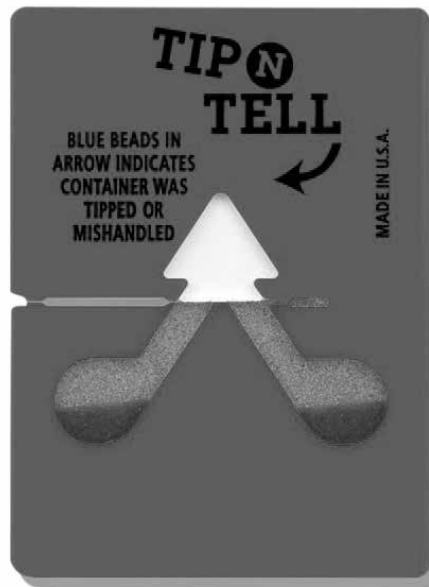

**Any damage should be reported immediately to the carrier and Peak Scientific or the Peak Partner from where the unit was purchased.**

Follow the unpacking instructions posted on the side of the crate. It will require two people to remove the unit from the shipping crate and to manoeuvre the generator to the desired location.

Please save the product packaging for storage or future shipment of the generator.

Note: Included with the generator is a "Fittings Kit" containing mains power leads for UK, EU & US and also all the required fittings and warranty registration card. Be careful not to discard these with the packaging.

## Fittings Kit Contents

Supplied in the Fittings Kit are all the fittings required to connect the generator to the application. The contents of the Fittings Kit are as follows:

1. Tube PE 6mmOD 4mmID	x3m
2. Cable C13 BS 1363 UK 230V 10A 2m	x1
3. Cable C13 Nema 6 15 US 220V 2m	x1
4. Cable C13 Schuko H05VV-F EU 230v 2.5m	x1
5. Cable C13 Nema 5 15 US 110V 2m	x1
6. Elbow 1/4"M BSPPx6mm Push Fit	x2
7. TOC 1000 Installation guide	x1

All of the generators output ports are located on the output panel at the rear of the unit.

# Installation

## Generator Environment

The generator should, if possible, be installed next to the application,. If, however, this is not possible, please consider the distance between the generator and your application, since pressure drops can result from extended lengths of tubing and may affect gas delivery.

Performance of the generator is affected by ambient conditions. Note should also be taken to the proximity of Air Conditioning outlets. Air conditioning systems can sometimes give rise to “pockets” of air with high relative humidity. Operation of the unit within such a pocket could adversely affect its performance. Consideration should also be given to the air flow around the unit. It is recommended that an air gap of 75mm (3”) should be maintained down both sides and at the rear of the unit. Please refer to the drawing below for the general dimensions of the unit.

Please ensure Generator is situated in a well ventilated environment and is positioned to allow easy disconnection from the mains supply if required.

Minimum Operating Ambient Temperature: 5 °C (41 °F)

Maximum Operating Ambient Temperature: 35 °C (95 °F)


**Compressed air supplied to the TOC 1000 generator should meet ISO 8573-1:2010 Class 1.4.1 standard and should be supplied at 100psi (6.9bar) - 145psi (10 bar) pressure. Supply of compressed air at a pressure of below 100psi may result in poor quality air that may not meet instrument requirements. Supply of compressed air at a pressure greater than 145psi may cause damage to the generator and could result in damage to any instrument supplied.**

**Compressed air supplied at <100psi could result in poor performance of the TOC 1000 gas generator and may affect instrument performance.**

**Compressed air supplied at >145psi may result in damage to the TOC 1000 gas generator and may result in damage to instrumentation connected to the gas generator.**


# Generator Overview

## General Dimensions

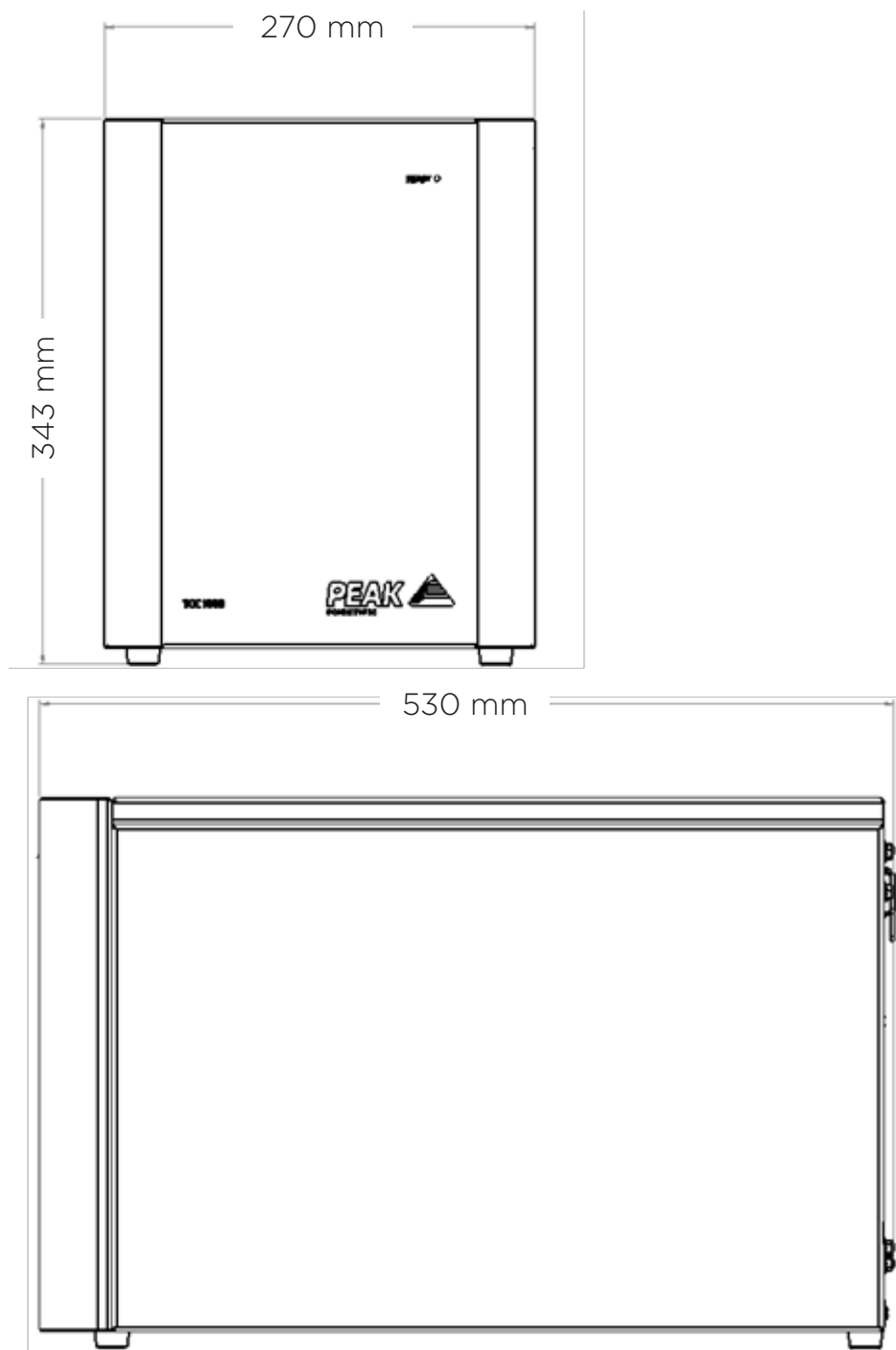

The generator must always be placed on a flat, level surface. Failure to do so will affect the performance of the generator.

## Rear Connections

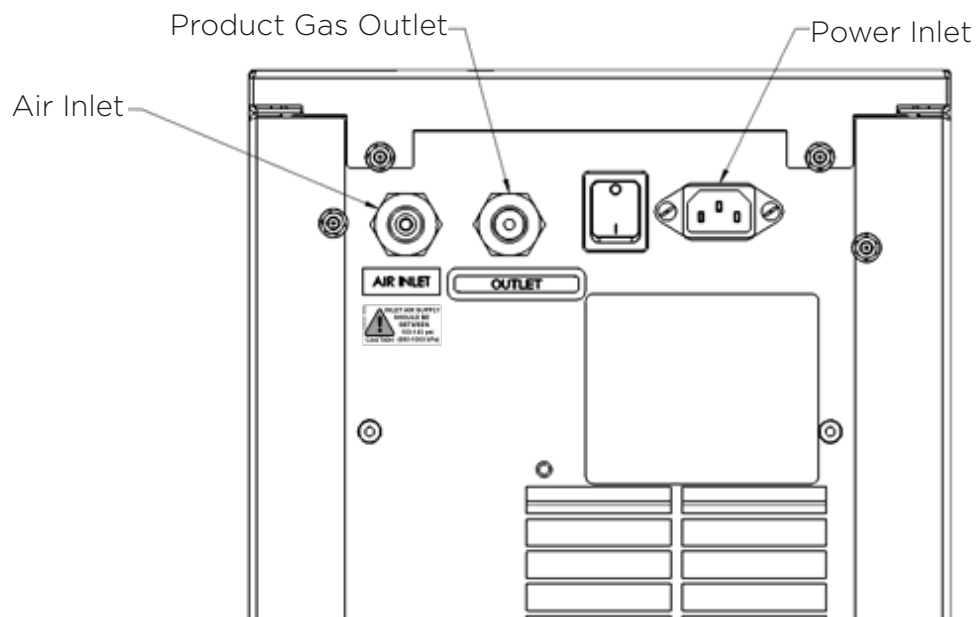

Ensure all inlets and outlets are connected to correct sources and applications.



All Connections should only be carried out by trained personnel. Generator must be switched off and unplugged prior to any cleaning or maintenance operations.



# Electrical Connection

Connect the generator to an appropriate 120 or 230 volt single-phase supply, refer to the generator serial plate for input specification and ensure your supply matches the requirements.

If the appropriate power cord is not supplied, a new plug with correct rating can be fitted by a qualified electrician, ensuring that it has adequate ratings and appropriate approvals for the country of operation. Failure to do so could cause damage to the generator or risk of overloading of the power cord.



If a substitute main supply cord is used, ensure that it has adequate ratings and appropriate approvals for the country of operation. Failure to do so could cause damage to the generator or risk of overloading of the power cord.

This unit is classified as SAFETY CLASS 1. THIS UNIT MUST BE EARTHED. Before connecting the unit to the mains supply, please check the information on the serial plate. The mains supply must be of the stated AC voltage and frequency.

EARTH/GROUND (E):-                      Green & Yellow                      or                      Green

LIVE (L):-                                      Brown                                      or                      Black

Neutral (N):-                                      Blue                                      or                      White

Electrical requirements are 120VAC nominal  $\pm 5\%$  or 230VAC nominal  $\pm 10\%$  depending on chosen model. Extended periods at extremes can have a detrimental effect on the operation and life of the generator.



If the equipment is used in a manner not specified by the manufacturer, the protection provided by the equipment maybe impaired.

## **Air Purity**

The Precision Zero Air Series generator should be connected to an air supply that, as a minimum, meets ISO8753-1:2010 Class 1.4.1

### **Class 1 Particulate**

In each cubic metre of compressed air, the particulate count should not exceed 20,000 particles in the 0.1 - 0.5 micron size range, 400 particles in the 0.5 - 1 micron range and 10 particles in the 1 - 5 micron size range.

### **Class 4 Water**

A pressure dew point (PDP) of +3°C or better is required and no liquid water is allowed.

### **Class 1 Oil**

In each cubic metre of compressed air, not more than 0.01mg of oil is allowed. This is a total level for liquid oil, oil aerosol and oil vapour.

## Start-Up Sequence



Before the generator is first connected to the application, the generator should be operated in isolation (i.e. not connected to the application) for 120 minutes. This is to ensure that any impurities are purged from the system. Failure to do this may harm the application.

## Connecting to the application

Using the tubing provided, connect the outlet of the unit to the inlet on the application.

If you require more tubing than is supplied please refer to the Tubing Lengths section below.



Once the tubing is connected to the application, please ensure that it is thoroughly checked for leaks. Even the slightest leak in the gas supply between the generator and the application can lead to a reduction in efficiency and gas purity.

**Compressed air supplied to the TOC 1000 generator should meet ISO 8573-1:2010 Class 1.4.1 standard and should be supplied air at a pressure of 100psi (6.9bar) - 145psi (10bar). Supply of compressed air at a pressure of below 100psi may result in poor quality air that may not meet instrument requirements. Supply of compressed air at a pressure greater than 145psi may cause damage to the generator and could result in damage to any instrument supplied.**

**Compressed air supplied at <100psi could result in poor performance of the TOC 1000 gas generator and may affect instrument performance.**

**Compressed air supplied at >145psi may result in damage to the TOC 1000 gas generator and may result in damage to instrumentation connected to the gas generator.**

## Tubing Lengths



The diameter of the tubing which will be connected to the gas outlet is important and is determined by the length of tubing required. Failure to follow these recommendations could lead to pressure drops between generator and application.

< 3 meters: Use 6mm OD / 4mm ID or 1/4" OD / 3/16" ID tubing.

> 3-10 meters: Use 8mm OD / 6mm ID or 5/16" OD / 1/4" ID tubing.

Some tubing and fittings not supplied in the fittings kit.

> 10 - 40 meters: Use 10mm OD / 8mm ID or 3/8" OD / 5/16" ID tubing.

> 40 metres: Please contact Peak Scientific with the relevant distance and we will calculate the flow resistance and the tubing size required.

# Service Requirements

## Service Schedule







Purchase Interval	Component	Qty.	Visit
12 Months	TOC 1000 Annual Maintenance Kit	1	<a href="http://www.peakscientific.com/ordering">www.peakscientific.com/ordering</a>

**Note: Only Peak approved components should be used.**

## Peak Protected

With Peak Scientific you invest in not only a product but peace of mind. With a network of certified Peak engineers stationed throughout the globe, Peak's rapid response team are never far away and our commitment is to keep your generator running day in, day out, protecting your laboratory workflow.
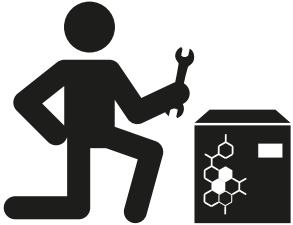
### [Peak Protected] can provide...

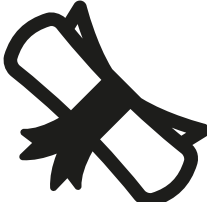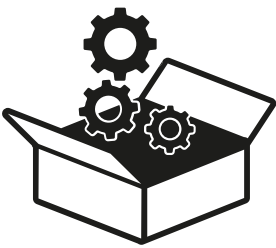 <p><b>Installation</b> A dedicated Peak engineer will visit your lab to install and setup your generator</p>	 <p><b>Complete plan</b> Swift response by a Peak Service Engineer within 72 hours &amp; planned preventative maintenance</p>	 <p><b>Premium Protected</b> Guaranteed rapid on-site response within 24 hours &amp; planned preventative maintenance</p>
 <p><b>IQ/OQ</b> Certified assurance for applications requiring documented qualification</p>	 <p><b>Spares</b> Genuine Peak parts with express delivery, ensuring optimal performance and lifetime</p>	 <p><b>Technical Support Hotline</b> Around the clock support by phone or online with our global technical helpdesk</p>

To find out more about protecting your investment visit: [www.peakscientific.com/protected](http://www.peakscientific.com/protected)

## Cleaning

Clean the outside of the generator only using warm soapy water and a clean damp cloth. Ensure all excess fluid is thoroughly removed from the cloth prior to use.



Cleaning should only be undertaken with the power switched off and the power cord removed from the rear of the generator.



Under no circumstances should any solvents or abrasive cleaning solutions be used as these can contain fumes that could be harmful to the generator.



Care should be taken with Leak Detection Liquids.

## Troubleshooting

Problem	Possible Solution
The application is reporting low pressure.	<ul style="list-style-type: none"><li>• Check for leaks between the TOC 1000 and the application.</li><li>• Contact your service provider.</li></ul>
Unstable baseline	<ul style="list-style-type: none"><li>• Check inlet air pressure is set to 100psi</li><li>• Check inlet air flow is 6LPM</li><li>• Check outlet air flow is 1LPM</li><li>• Contact your service provider</li></ul>



# Go Online or Complete and Return


We know that registering any of your recently purchased products is not the first thing on your mind- but it is very important to both of us. Not all warranties are alike and Peak Scientific stand out against other gas suppliers as we offer a comprehensive, quick response, on-site warranty. This means that in the very unlikely case that your gas generator develops a fault we have rapid support teams on-hand around the world who are able to come to your lab and get you back up and running in no time.

Register for your **comprehensive 12 month on-site warranty** with ease online at [www.peakscientific.com/protected](http://www.peakscientific.com/protected).

Alternatively, you can send the completed form to Peak Scientific by post or email at [warranty@peakscientific.com](mailto:warranty@peakscientific.com).



The screenshot shows the 'PEAK' logo at the top left. Below it, the 'Serial number' section includes a 'Generator Serial Number' field with a dropdown arrow and an 'Installation date' field with a calendar icon. The 'Contact details' section contains fields for 'First name', 'Last name', 'Email address', and 'Phone'. A 'Submit' button is at the bottom right.

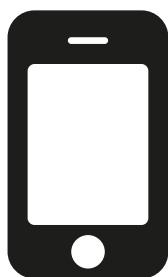

The screenshot shows a mobile interface with a green header and the title 'Your Details'. It includes input fields for 'First name', 'Last name', 'Email address', 'Phone number', and 'Address'. A 'City' field is at the bottom. A 'Submit' button is in the bottom right corner.

## Important!

You have **1 month to register** your Peak Scientific product from the date of installation. Once registered the warranty will be honoured for a period of 12 months. If you wish to defer the installation of your generator, you must notify Peak Scientific immediately by emailing **[warranty@peakscientific.com](mailto:warranty@peakscientific.com)**. For generators that remain unregistered after 1 month from the shipment date, the warranty will be considered active from the date of factory dispatch.





# [**PEAK Protected**]<sup>TM</sup>

Peak Scientific has highly trained, fully certified Field Service Engineers located in over 20 countries across every continent around the world. This allows us to provide an industry-leading rapid response service to our customers. With **[Peak Protected]**, your laboratory's productivity becomes our top priority.

To discuss Peak Protected generator cover and payment options speak to your local Peak Representative or for further information contact: **[protected@peakscientific.com](mailto:protected@peakscientific.com)**

## **Peak Scientific**

Fountain Crescent  
Inchinnan Business Park  
Inchinnan  
PA4 9RE  
Scotland, UK

**Tel:** +44 141 812 8100

**Fax:** +44 141 812 8200

For further information on any of our generator products please contact **[discover@peakscientific.com](mailto:discover@peakscientific.com)**

