

Your local *gas generation* partner



Gas Generators

For Agilent Instruments



Agilent Technologies

Discover more at www.peakscientific.com/agilent

Our story

PEAK PERFORMERS: OUR WORLD-CLASS PRODUCT PLATFORMS



Four distinct product series define our core offering to supply Agilent Technologies analytical laboratory instruments. Our bestselling Genius and Genius XE Series, designed mainly for LC-MS/MS applications, is a self-contained system featuring integrated compressors, delivering laboratory-grade nitrogen at various pressures and flow rates. We also have a specific solution for Agilent MP-AES instruments, the Halo air & nitrogen generator launched in 2018.

The Infinity, Infinity XE 60 & Solaris XE Series of nitrogen gas generators are silent-running compressor-free systems for laboratories with an in-house supply of clean, dry air – delivering nitrogen in a wide range of flow rates.

Engineered specifically to meet the needs of Gas Chromatography our cutting-edge Precision Series is a modular system consisting of high-purity nitrogen, hydrogen and zero air generators that can be combined into a compact stack tailored to your needs. New Precision Series models featuring custom blue LED lighting have been introduced exclusively for Agilent.

For very high-flow, high demand needs, our innovative Pressure Swing Adsorption based i-FlowLab technology is designed to meet the full and total demand for nitrogen gas a facility may have, all in one expandable installed solution.

OUR PROMISE

Our commitment is to provide laboratories with a reliable gas generation solution that delivers long term value. We achieve this through exceptional product design, and direct global on-site support that is unmatched in our industry.

We take a vested interest in customer outcomes and strive to be a valued partner in their success by delivering exceptional service over the life of our product.

Our story

THE PEAK WAY



PEAK Scientific is a leading innovator in the design, manufacture and support of high-performance gas generators for analytical laboratories. Established in 1997 near Glasgow (Scotland), where our corporate Headquarters and high-tech manufacturing and R&D facilities reside, PEAK Scientific boasts a significant local presence on every continent – including major operations in North America, China and India.

With over two decades of experience pioneering reliable gas generator technology, PEAK Scientific develops market-leading nitrogen, hydrogen and zero air systems mainly for the field of LC-MS and GC.

Our strength lies in providing our customers with high-quality, dependable gas generation solutions tailored to meet their specific laboratory demands; what differentiates us is our world-class technical support and on-going service care throughout the generator's lifespan, wherever you are in the world.

MANUFACTURING EXCELLENCE



Our products are the result of our meticulous Research and Development culture, built over more than 20 years at the forefront of the gas generator market. We pride ourselves on the greatest care taken to assess specific application needs before designing and rigorously testing new products.

As we have expanded, our R&D capabilities have advanced to the point where we are better equipped than ever to meet changing market demands. By implementing the latest manufacturing technologies and philosophies, PEAK continues to improve product quality, responsiveness and efficiency – resulting in a better value proposition for our customers.

Every PEAK Scientific generator is designed from our ISO 9001:2015 accredited manufacturing center of excellence in the UK and independently tested to ensure compliance with all applicable safety standards.

Application Matrix

| | | PEAK Solution | |
|--------------------------------------------------|---------------------|----------------------------------------|----------------------------------------|
| Application | | Single Instrument | Two Instruments |
| 1200 Series Infinity ELSD | No compressor | Solaris 10 | Solaris 10 |
| | Built in compressor | Solaris 10 with Solaris Air Compressor | Solaris 10 with Solaris Air Compressor |
| 385 ELSD | No compressor | Solaris 10 | Solaris 10 |
| | Built in compressor | Solaris 10 with Solaris Air Compressor | Solaris 10 with Solaris Air Compressor |
| InfinityLab MSD-iQ | No compressor | Solaris XE 35 | Solaris XE 35 |
| | Built in compressor | • Genius SQ 24 | • Genius XE 70 |
| LC-MS with ESI | No compressor | Solaris XE 35 | Solaris XE 35 |
| | Built in compressor | • Genius SQ 24 | • Genius XE 70 |
| LC-MS with APCI | No compressor | Solaris XE 35 | Solaris XE 35 |
| | Built in compressor | • Genius XE 35 | • Genius XE 70 |
| All LC-MS systems / sources except 6490/95, 6550 | No compressor | Solaris XE 35 | Infinity XE 6010 |
| | Built in compressor | • Genius XE 35 | • Genius XE 70 |
| LC-MS 6490/95, 6550 | No compressor | Infinity XE 6010/1 | - |
| | Built in compressor | • Genius XE 70 | - |
| ICP-OES 700 Series | No Compressor | • Precision Nitrogen NG Range | • Precision Nitrogen NG Range |
| ICP-MS 7500cs/ 7500cx | No Compressor | • Precision Hydrogen Trace 250cc | • Precision Hydrogen Trace 250cc |
| ICP-MS 7700x | No Compressor | • Precision Hydrogen Trace 250cc | • Precision Hydrogen Trace 250cc |
| ICP-MS 7900 | No Compressor | • Precision Hydrogen Trace 250cc | • Precision Hydrogen Trace 250cc |
| Cary Eclipse Fluorescence Spectrophotometer | No Compressor | i-FlowLab 701X | i-FlowLab 701X |
| Cary 600 Series FTIR Spectrometer | No Compressor | PG14L Precision compressor | PG28L |
| Cary 610/20 Series FTIR Microscope | No Compressor | PG14L | PG28L |
| • Included on Agilent Corporate Price List | | | |

PEAK Gas Generator Recommendations For Agilent GC

| Detector | Carrier Gas | | Detector Gas | | |
|-----------------|----------------------------|--------------------------|--------------------------|--------------------|--------------------------|
| | Hydrogen | Nitrogen | Hydrogen | Air | Nitrogen |
| MS | • Precision Hydrogen Trace | Precision Nitrogen Trace | Not Applicable | | |
| FID | • Precision Hydrogen Trace | Precision Nitrogen Trace | Precision Hydrogen | Precision Zero Air | Precision Nitrogen |
| TCD | • Precision Hydrogen Trace | Precision Nitrogen Trace | - | - | Precision Nitrogen |
| ECD | • Precision Hydrogen Trace | Precision Nitrogen Trace | - | - | Precision Nitrogen Trace |
| ECD Sensitivity | • Precision Hydrogen Trace | Precision Nitrogen Trace | - | - | Precision Nitrogen Trace |
| FPD | • Precision Hydrogen Trace | Precision Nitrogen Trace | Precision Hydrogen Trace | Precision Zero Air | - |
| FTD | • Precision Hydrogen Trace | Precision Nitrogen Trace | Precision Hydrogen | Precision Zero Air | Precision Nitrogen |

• Included on Agilent Corporate Price List



Product

Genius XE Series

DESCRIPTION

Inspired by the success of our best-selling Genius line of nitrogen gas generators for LC-MS-MS, Genius XE Nitrogen is a cutting-edge evolution combining advanced technology with refined and robust engineering. With two models - XE 35 (up to 35 L/min) and XE 70 (up to 70 L/min) - Genius XE Nitrogen provides a premium standalone nitrogen solution for high performance LC-MS-MS and other mission-critical laboratory applications where performance and reliability are paramount.

Features

Variable flow up to 70 L/min

Variable pressure up to 116 psi

Multi-Stage Purification™ producing analytical grade nitrogen gas up to 99.5% purity

2 year comprehensive manufacturer's warranty*

Touch-screen full colour user interface for ease of operation

Additional compressor capacity (Genius XE 70 only)



| Product | Flow Rate | Gas Output | Pressure | Size (HxWxD) | Application |
|--------------|-----------|------------|--------------|--------------------------------------------|------------------|
| Genius XE 35 | 35 L/min | Nitrogen | 116psi /8bar | 650 x 570 x 710 mm 25.2 x 22.4 x 28.0" | LC-MS-MS or ELSD |
| Genius XE 70 | 70 L/min | Nitrogen | 116psi /8bar | 1000 x 570 x 710 mm 39.4 x 22.4 x 28.0" | LC-MS-MS or ELSD |

* Year 2 of warranty subject to generator being serviced at end of year 1 by a PEAK-approved agent in accordance with fixed annual maintenance schedule.



Product

Genius SQ 24

DESCRIPTION

For labs using any Single Quad LC-MS systems, the Genius SQ 24 nitrogen gas generator delivers an analytical grade nitrogen supply dedicated for this instrument. With a compact size to fit under most lab benches and flow rates of up to 24 L/min, meeting and exceeding the maximum flow rate required by any Single Quad LC-MS on the market today, your lab can enjoy hassle free nitrogen gas which helps you deliver the best analytical results.

Features

Best value nitrogen generator on the market

Variable flow up to 24 l/min

Analytical grade nitrogen supply for Single Quad LC-MS

Outlet pressure up to 116psi (at 22 l/min maximum)

Reduced height and compact size to fit under most lab benches

Integrated internal compressors for easy plug and play



| Product | Flow Rate | Gas Output | Pressure | Size (HxWxD) | Application |
|--------------|-----------|------------|--------------|-------------------------------------------|-------------------|
| Genius SQ 24 | 24 L/min | Nitrogen | 116psi /8bar | 610 x 600 x 750 mm 24.0 x 23.6 x 29.5" | Single Quad LC-MS |

Product

Solaris XE 35

DESCRIPTION

PEAK Scientific's Solaris XE 35 generator has been engineered using membrane technology to meet the gas delivery requirements for labs using LC-MS or multiple lower flow nitrogen instruments such as ELSDs or Compact Mass Spectrometers. Solaris XE can operate with varying flow rates (up to 35L/min of high purity nitrogen), purity (up to 99.5%) and outlet pressure can be adjusted down from 116psi, offering a flexible solution for a variety of applications.

Features

Nitrogen up to 35L/min, up to 116 psi and up to 99.5% purity

Variable outlet flow, pressure and purity

Compact, stackable system maximising the use of valuable laboratory floor or bench space

Colour indicative LED lighting for easy status awareness

Manufactured and tested to highest spec in UK, CE/FCC/CSA certified



| Product | Flow Rate | Gas Output | Pressure | Size (HxWxD) | Application |
|---------------|-----------|------------|--------------|------------------------------------------|-----------------------|
| Solaris XE 35 | 35 L/min | Nitrogen | 116psi /8bar | 156 x 343 x 650 mm 6.1 x 13.5 x 25.6" | LC-MS, Multiple ELSDs |

Product

Infinity XE 6010

DESCRIPTION

Producing between 10 and 500 L/min of high purity nitrogen gas, the Infinity XE 60 series nitrogen membrane gas generators can comfortably supply multiple laboratory instruments with instrument grade nitrogen.

Features

- Flow rates ranging from 10-500 L/min*
- N₂ purity up to 99.5%
- Versatile solution with varying Nitrogen flow, pressure & purity capabilities
- In-field expansion of nitrogen capacity to meet growing lab demands
- Optional second outlet for Dry Air
- Intelligent Membrane Management (IMM) optimizes compressed air requirement based on gas consumption



| Product | Flow Rate | Gas Output | Pressure | Size (HxWxD) | Application |
|------------------|-----------------------------------|---------------------|------------------------------------|-------------------------------------------|--------------------------------------|
| Infinity XE 6010 | up to 500 L/min up to 210L/min | Nitrogen Dry Air | 135 psi / 8 bar 145psi / 10 bar | 995 x 410 x 950 mm 39.2 x 16.1 x 37.4" | For multiple LC-MS/MS instruments |

Product

Precision Hydrogen

DESCRIPTION

The Precision Hydrogen generators come in two variants, standard and trace. The standard models are designed to provide detector gas for applications such as FID available in 100cc, 200cc, 300cc, 450cc and 1200cc/min. The trace model can deliver 250cc, 500cc or 1200cc/min of UHP hydrogen for both carrier and detector gas. One generator is capable of supplying multiple detectors and GC or GC-MS units.

Features

99.9995% purity (standard) / 99.99999%* purity (trace model)

Suitable for GC detector gas (all) and carrier gas (trace model only)

Creates hydrogen on demand, minimal storage of hydrogen in the system

Low maintenance and minimal running costs over product lifetime

Internal leak detection with automatic shutdown

Automatic loading pump as standard

** Based on O₂ content independently verified by National Physical Laboratory, UK*



Product

Precision Nitrogen

DESCRIPTION

The Precision Nitrogen generators are available in both standard and trace models. The standard models come in 250cc, 600cc and 1 liter and are able to provide make-up gas, reference gas and for sample preparation (i.e. headspace). The Nitrogen Trace models are engineered to offer zero nitrogen as carrier gas, make-up gas and reference gas for sample preparation available in 250cc, 600cc and 1L models.

Features

All models offer 99.9995% purity

Nitrogen Trace model suitable as GC and GC-MS carrier and detector gas

Trace capable of delivering hydrocarbon-free Nitrogen (as Methane) suitable for carrier gas, make-up gas and sample preparation

Highly economical source of nitrogen gas with low lifetime running costs

Ultra-fast start-up time, quick to reach standard operating purity



Product

Precision Zero Air

DESCRIPTION

The Precision Zero Air generators are designed specifically to supply clean, dry, hydrocarbon-free air to be used as flame support gas for GC.

As with all Precision Series generators, Zero Air benefits from a compact and modular, stackable design, minimizing the total footprint required for GC gas supply, and providing flexibility to add or remove modules as your laboratory requirements evolve over time.

Features

Hydrocarbon content <0.05ppm for market-leading purity

Integrates seamlessly with other Precision units

Minimum lifetime maintenance requirements, no expensive catalyst chamber replacements

Highly economical source of dry, hydrocarbon free air

Avoid risk of contaminants entering the system (when switching out empty cylinders)



Peak Protected

Safeguard your laboratory workflow with [**PEAK** Protected]TM

Servicing any equipment in your lab should not be seen as just an option, it should be treated as a necessity. With many moving parts and wear & tear, especially where generators feature integrated compressors, as well as filters and other 'consumables', a gas generator needs to be maintained regularly to continue to perform at its best for many years.

Here at PEAK we see your gas generator as the beating heart of your laboratory, enabling your analytical instruments to perform and achieve the results you need. With this being the case it is extremely important you look after your 'beating heart' to ensure the best possible results.

We offer a range of **[PEAK Protected]** services to suit your particular needs.



Premium Plan

Complete maintenance care with 24 hour rapid response breakdown cover



Complete Plan

Complete maintenance care with guaranteed response time breakdown cover



Fixed Price PM

Preventative maintenance carried out by an expert PEAK engineer at a time that suits you



Supported Plan

Preventative maintenance support specifically for compressor-less generators



Paid Repairs

On-site repairs carried out by a certified PEAK FSE in case of a breakdown



Replacement Parts

Genuine PEAK parts with express delivery, ensuring optimal performance and lifetime



Installation

Trained PEAK FSE will visit to install and commission your generator



IQ/OQ

Certified assurance for applications requiring documented qualification



Support Hotline

Around the clock support by phone or online with our global technical helpdesk


Not all warranties are equal.

What differentiates PEAK from other gas generator manufacturers is that a PEAK gas generator comes with a truly comprehensive on-site warranty. This means that in the unlikely event that your gas generator develops a problem, we will send a PEAK Certified Field Service Engineer to your laboratory to try to resolve the issue on-site and get you back up and running with minimal fuss and disruption.

What's more, when we say fully-comprehensive, we mean it. Every part of your generator is covered during warranty. While other manufacturer warranties may only cover certain components or require the return of your generator to the factory for repair, leaving you without gas for weeks or even months, our warranty is designed to give you complete peace of mind.

To find out more about our unique warranty, visit www.peakscientific.com/warranty

Your local **gas generation** partner

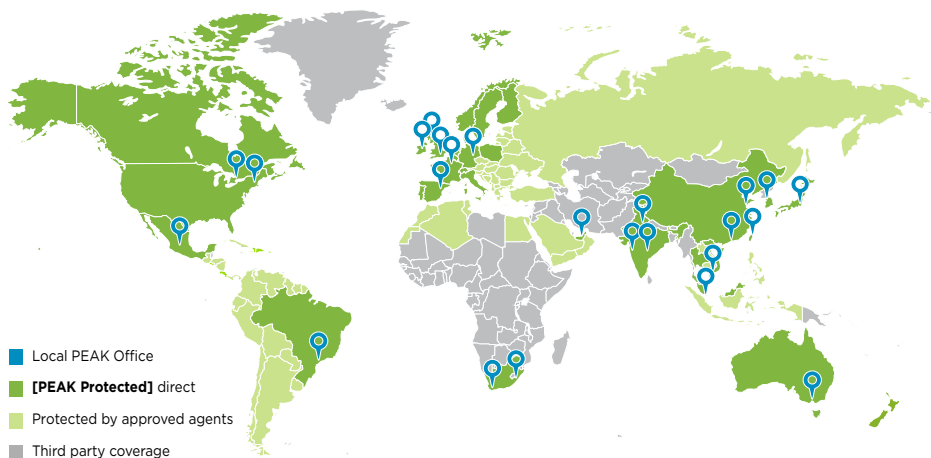

PEAK Scientific has highly trained, fully certified Field Service Engineers located in over 20 countries across every continent around the world. This allows us to provide an industry-leading rapid response service to our customers. With **[PEAK Protected]**, your laboratory's productivity becomes our top priority.

Contact us today to discover more!

North America
Tel: +1 866 647 1649

Europe
Tel: +44 (0)141 812 8100

China
Tel: +86 21 5079 1190

India
Tel: 1800 2700 946

For contact details of all PEAK local offices around the world, visit:
www.peakscientific.com/offices